
CHOLESTEROL OXIDASE [CON II]

Liquid type

from *Rhodococcus* sp.

(Cholesterol: oxygen oxidoreductase, EC 1.1.3.6)

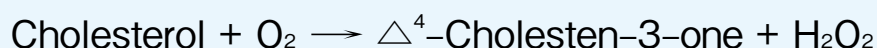Preparation and Specification

Appearance : Yellowish solution
Specific activity : More than 500 U/ml

Properties

Substrate specificity	: See Table 1	
Molecular weight	: 61.8 KDa (SDS-PAGE)	
Isoelectric point	: pH 4.5	
Michaelis constant	: Cholesterol $6.0 \times 10^{-5}\text{M}$	
Optimum pH	: 7.0-7.5	Figure 1
pH stability	: 5.7-7.8 (65°C, 10 min)	Figure 2
Optimum temperature	: 50°C (Tris-HCl buffer)	Figure 3
Thermal stability	: Stable at 65°C and below (pH 7.0, 10 min)	Figure 4 and Figure 5
Effects of detergents	: See Table 2	

Applications for Diagnostic Test

This enzyme is useful for enzymatic determination of **total cholesterol**, **HDL-C**, and **LDL-C** when coupled with cholesterol esterase (T-18 and T-98).


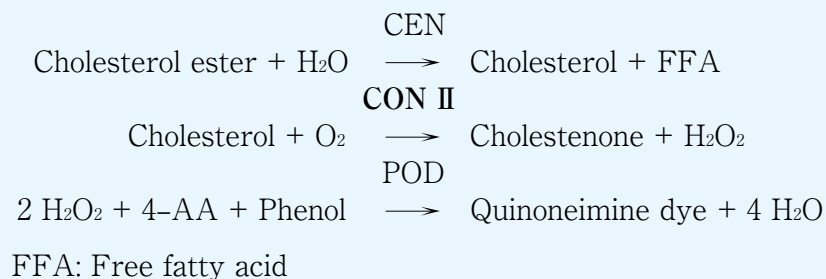
Table 1. Substrate specificity

Substrate (1mM)	Relative activity (%)
Cholesterol	100
β -Cholesterol	93
Pregnenolone	98
Dehydro-iso-androsterone	10
β -Sitosterol	94
Stigmasterol	66
Androsterol	2
Teststerone	1
Cholic acid	3

Table 2. Effect of detergents on CON II activity


Detergents (0.1%)	Relative activity (%)
Triton X-100	100
Emulgen 810	101
Emulgen 911	113
Emulgen 709	107
Emulgen 109P	118
AdekatoI B-797	103
AdekatoI SO-120	100
AdekatoI 720N	114
RHEODOL 460	63
SM 1080	122

Fig.1 pH Optimum

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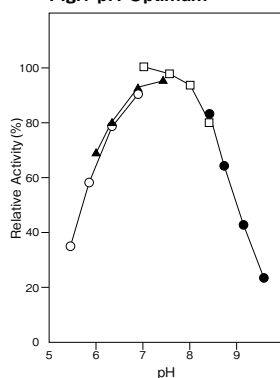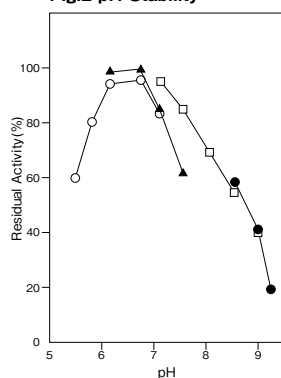
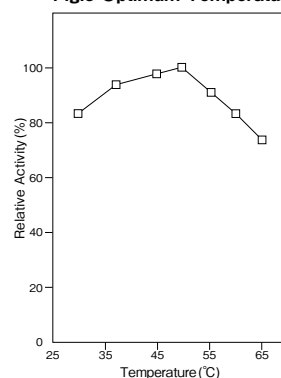
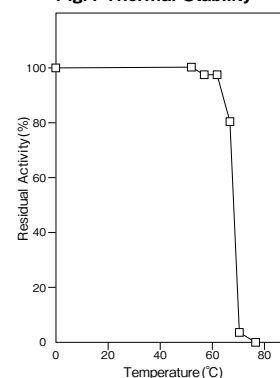
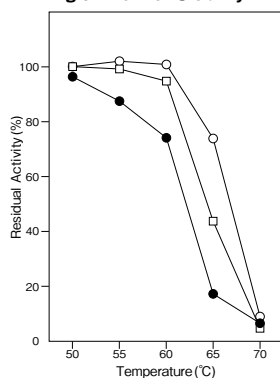
Assay

Principle

The assay is based on the increase in absorbance at 240 nm as Δ^4 -cholesten-3-one is produced in the following reactions:

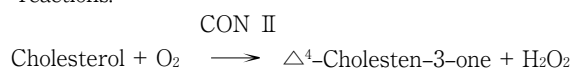Unit definition

One unit is defined as the amount of enzyme which liberates 1 μ mole of Δ^4 -cholesten-3-one per minute at 37°C under the conditions specified in the assay procedure.

Reagents

- Substrate solution (6 mM cholesterol solution)
 Dissolve 232 mg of cholesterol with isopropanol to make a total of 100 ml.

2. Enzyme dilution buffer
0.1 M KH_2PO_4 - Na_2HPO_4 buffer pH 7.0 containing
0.05% (W/V) Triton X-100
3. Reagents
Cholesterol : NACALAI TESQUE, INC. Special grade
#08721
Triton X-100 : The Dow Chemical Company

■ Enzyme solution

Dilute accurately 0.5 ml of the sample with enzyme dilution buffer to make a 40-fold solution. Dilute it with enzyme dilution buffer to adjust the concentration to within 0.1-0.2 U/ml.

■ Procedure

1. Pipette accurately 3.0ml of enzyme dilution buffer and 50 μl of enzyme solution and preincubate at 37°C.
※ In the case of a test blank, add 50 μl of enzyme dilution buffer in place of enzyme solution.
2. After 5 min, add 50 μl of substrate solution and mix to start the reaction at 37°C.
3. After starting the reaction, measure the rate of increase per minute in absorbance at 240nm. The rate must be measured within the linear portion of the absorbance curve.

Absorbance sample : A_s/min
blank : A_b/min

$$0.010 \text{ Abs/min} \leq \Delta A/\text{min} = (A_s/\text{min} - A_b/\text{min}) \leq 0.060 \text{ Abs/min}$$

■ Calculation

$$\text{Activity (U/ml)} = \frac{\Delta A/\text{min}}{12.2} \times \frac{3.10}{0.05} \times D$$

- 12.2 : millimolar extinction coefficient of Δ^4 -Cholesterol-3-one at 240 nm ($\text{cm}^2 / \mu\text{mole}$)
3.10 : final volume (ml)
0.05 : volume of enzyme solution (ml)
D : times of dilution in enzyme solution

Storage

Storage at -80°C in the presence of a desiccant is recommended. Enzyme activity will be retained for at least one year under this condition.

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