



POINT- OF-CARE

osom[®]

Ultra hCG Combo Test

OSOM[®] Ultra hCG Combo Test

FOR THE DETECTION OF HUMAN CHORIONIC GONADOTROPIN IN URINE OR SERUM

Features and Benefits

- **Accurate** – No biotin interference. Detects beta-core fragments and other HCG variants*.
- **Fast and Easy** – Single step procedure with results in 3 minutes
- **Convenient** – Pipette included in each test pouch
- **Exceptional Customer Service** – Online training modules and technical support staffed by certified medical technologists/professionals

Performance Characteristics

Sample Type	Limit of Detection	Storage Temperature	Time to Result
Serum, urine**	25 mIU/mL	2-30°C/36-86°F	3 minutes


* Data on file

**CLIA Waived for urine



OSOM® Ultra hCG Combo Test

SIMPLE TO USE / RESULTS IN MINUTES

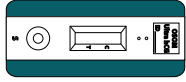


1 ADD 3 DROPS OF SAMPLE INTO THE SAMPLE WELL



2 READ THE RESULTS AT 3-5 MINUTES

POSITIVE

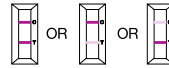

NEGATIVE



INVALID



Ordering Information

	Configuration	Catalog Number	Storage	CPT Code
OSOM ULTRA HCG COMBO TEST	25 Test + 2 test for QC	1004	Room Temp	Urine: 81025QW Serum: 84703
OSOM HCG URINE CONTROL SET	2 X 10mL Vials (Pos & Neg)	134	Refrigerated	
OSOM HCG SERUM CONTROL SET	2 X 5mL Vials (Pos & Neg)	138	Refrigerated	



To learn more about all of our OSOM Tests scan the code or go to [osomtests.com](https://www.osomtests.com)



THE AMERICAS

SEKISUI Diagnostics, LLC
One Wall Street
Burlington, MA 01803
Phone: 800 332 1042
Fax: 800 762 6311
Email: questions@sekisui-dx.com

SEKISUI
DIAGNOSTICS

Because every result matters™

[sekisuidiagnostics.com](https://www.sekisuidiagnostics.com)

